

SIEMENS



www.siemens.com/verio

Upgrade your MAGNETOM Verio
to the new MAGNETOM Skyra^{fit}

Transforming 3T productivity.

Answers for life.

Upgrade now

and benefit from the latest innovations in 3T

You are working with MAGNETOM Verio – a perfect starting point to take the next step in MRI innovation. A powerful upgrade is now available for your system, making your 3T imaging more accurate, more productive, and more consistent. This will help you to stay secure and ahead with your investment. Today, and for years to come. Making sure patients benefit most in the end.

Your proven MAGNETOM® Verio can now be upgraded with the key technologies Tim® 4G and Dot®. The new MAGNETOM Skyra^{fit} is now fit for the future, fitting your changing needs in a dynamic healthcare environment. Your 3T magnet will be working with the 45 mT/m @ 200 mT/s VQ Gradients together with Tim 4G, the 4th generation of proven Tim (Total imaging matrix) integrated coil technology. Tim 4G – with its all new RF architecture – locates all transmit and receive components at the magnet, resulting in a fully digital architecture. Your benefits include exceptional image quality, faster acquisition and exam time, as well as higher signal-to-noise ratio (SNR).

The upgrade includes Dot (Day optimizing throughput), the next movement in MRI. Dot is a new way of scanning in MRI – a better way. Dot scanning uses a suite of customizable engines – allowing the user to personalize exams according to patient needs, build in step-by-step user guidance, and automate MRI exams – either “out of the box” or based on the institution’s standards. Your benefits include increased consistency and reproducibility, greater ease-of-use, and higher productivity.



Your MAGNETOM Verio		The new MAGNETOM Skyra ^{fit}
MAGNETOM Verio Magnet	>>	MAGNETOM Verio Magnet remains
VQ Gradients 45 mT/m @ 200 T/m/s	>>	VQ Gradients 45 mT/m @ 200 T/m/s
Tim [102x8], [102x18], [102x32]	>>	Tim 4G [204x48], [204x64]
syngo MR B17	>>	syngo MR D13 with new advanced applications
Without Dot	>>	Up to 9 Dot engines

Upgradeability. With MAGNETOM.

Upgradeability solutions for each business need. The new system upgrade including Tim 4G and Dot will bring your system to the latest innovative hardware and software platform in MRI.



Step-by-step from MAGNETOM Verio to MAGNETOM Skyra^{fit}

Based on your original 3T magnet, you will immediately benefit from MAGNETOM Skyra^{fit} latest technologies: Exciting new applications, Siemens' revolutionary Tim 4G architecture, Dual-Density Signal Transfer,

and the Dot MRI scanning software platform. And: The upgrade of your MAGNETOM Verio to MAGNETOM Skyra^{fit} can be completed in only up to 15 working days!



Magnet room

The Body Coil is removed and replaced with a new one.



Licenses migration

Installed licenses are migrated into syngo[®] MR D.¹



¹ Consult your local Siemens representative for further details.




Technical room

Control unit and cooling unit cabinets are removed and replaced with new ones.



Operator's room

All workstations, monitors, and keyboards are removed and replaced by new ones.

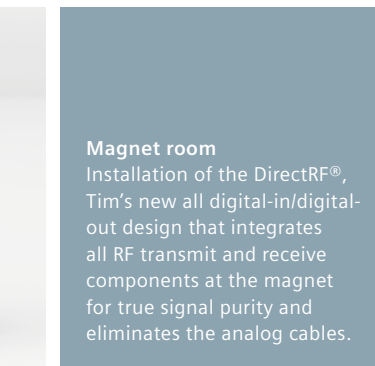

Magnet room

Installation of the DirectRF®, Tim's new all digital-in/digital-out design that integrates all RF transmit and receive components at the magnet for true signal purity and eliminates the analog cables.



Magnet room

All removed covers are replaced by new ones with Dot Display and Dot Control Centers. The patient table is replaced with the new Tim Table or the optional Tim Dockable Table.




Hand-over

After installation and image quality test, a comprehensive application training is held to help you get the best out of your new system.



More possibilities
for your 3T magnet



Up to 50% more
patient throughput
with Tim and Dot.¹



Are you ready to take the
next step in 3T innovation?

Time to

¹ Data on file. Results may vary.

change



“Dot has the advantages of an automatic gear box: Changing gears is unnecessary, but we stay at the wheel and decide where to drive.”

*Dr. Professor Henrik Michaelis, Section Chief
Vascular and Abdominal Radiology, Institute for
Clinical Radiology and Nuclear Medicine, Uni-
versity Medical Center Mannheim, Germany*

This material is preliminary and for internal use only. The content has not undergone regulatory review.

Higher reproducibility, higher productivity

Dot, the next movement in MRI, is a new way of scanning in MRI. Dot scanning uses a suite of customizable engines – allowing the user to personalize exams according to patient needs, build in step-by-step user guidance, and automate MRI exams. Your benefits include increased consistency and reproducibility, greater ease-of-use, and higher productivity. With MAGNETOM Skyra^{fit}, the following Dot engines are available:

Brain Dot Engine

- Up to 20% higher throughput and faster reading¹
- Reproducible positioning and standardized image quality
- Consistent exam duration and more efficient scheduling



Abdomen Dot Engine

- Up to 28% better timing accuracy¹
- Consistent image quality for even complex abdomen examinations



Cardiac Dot Engine

- Up to 50% increase in patient throughput¹
- Consistency in slice positioning for reliable image quality over multiple exams
- Ease-of-use helps bringing cardiac MRI into your clinical routine



¹ Data on file. Results may vary. ² MR imaging of patients with metallic implants brings specific risks. However, certain implants are approved by the governing regulatory bodies to be MR conditionally safe. For such implants, the previously mentioned warning may not be applicable. Please contact the implant manufacturer for the specific conditional information. The conditions for MR safety are the responsibility of the implant manufacturer, not of Siemens.

Angio Dot Engine

- Automated calculation of contrast agent application
- Interactive contrast timing approach eliminates need for cumbersome calculations
- Increased timing accuracy and image consistency



TimCT Angio Dot Engine

- Combines integrated Continuous Table move with Dot's optimized timing
- Up to 30% faster acquisitions for whole-body angio exams¹



TimCT Onco Dot Engine

- Combines integrated Continuous Table move with Dot's dynamic acquisition of upper abdomen
- Scan thorax, abdomen, and pelvis in less than 30 minutes¹



Breast Dot Engine

- Faster planning
- Higher standardization
- Less errors



Spine Dot Engine

- Complete spine examinations with ease
- Fast and standardized scanning
- Consistent and robust image quality




Large Joint Dot Engine

- Standardized examination and reproducible positioning with AutoCoverage and AutoAlign
- syngo WARP for reduction of susceptibility artifacts, such as from MR conditional metal implants²




New possibilities, new applications

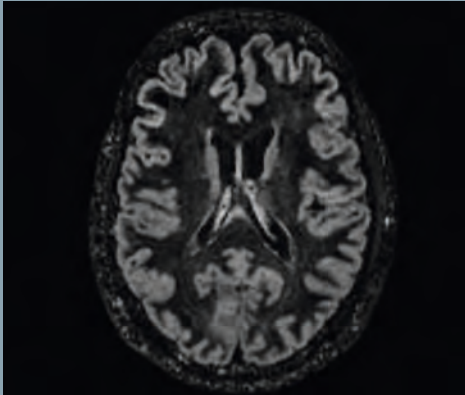


***syngo* SPACE DIR**

New with *syngo*® MR D13 are Double Inversion Recovery 3D protocols (SPACE DIR) with two user selectable inversion pulses for the simultaneous suppression of e.g. cerebro-spinal fluid and white matter.

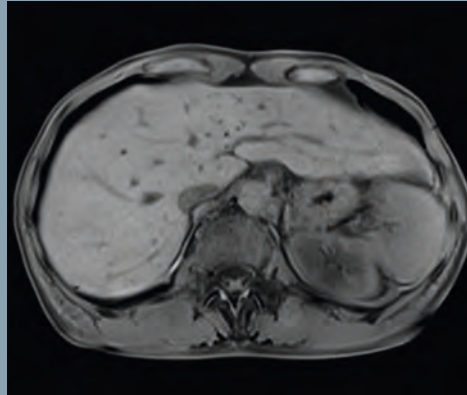

CAIPIRINHA


A new iPAT² sequence technique named CAIPIRINHA (Controlled Aliasing In Parallel Imaging Results IN Higher Acceleration) has been added. It can be applied to volumetric 3D imaging e.g. in the abdominal region. Higher PAT factors require more oversampling. CAIPIRINHA pattern distributes k-space points more uniformly.

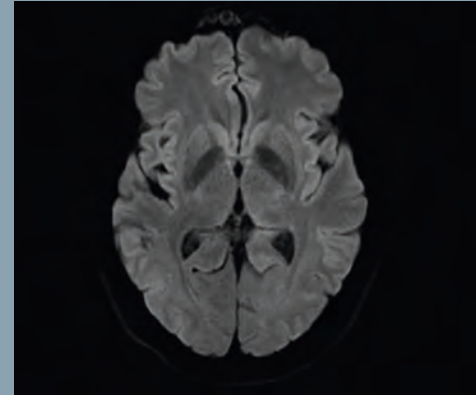

***syngo* RESOLVE**

syngo RESOLVE (Readout Segmentation Of Long Variable Echo-trains) delivers high-resolution Diffusion-Weighted Imaging (DWI) to visualize the diffusion properties of fine anatomical structures, enabling accurate lesion evaluation.

Additionally, this technique is largely insensitive to susceptibility effects, providing detailed anatomy-true diffusion imaging for brain, spine, breast, and prostate.

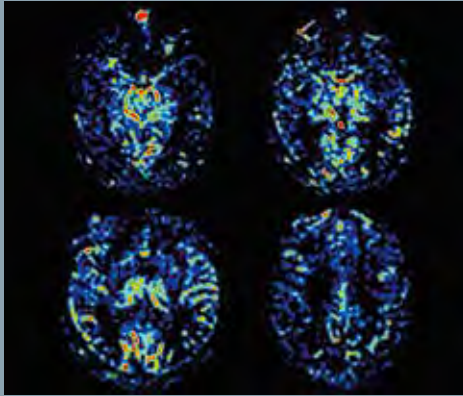


syngo ASL 3D

3D Arterial Spin Labeling (ASL) allows the non-invasive evaluation of brain perfusion without the injection of a contrast agent.


syngo MR D13 offers the possibility of measuring 3D ASL with multiple inversion times (TI).

This option provides data which can be used for Bolus Arrival Time (BAT) map calculations. Multiple TI acquisition allows for acquisition of raw label-control pairs at different (equidistant) inversion times. 3D ASL is the radiation-free alternative to PET as it offers an increase in SNR and a shorter scan time with reduced motion sensitivity.

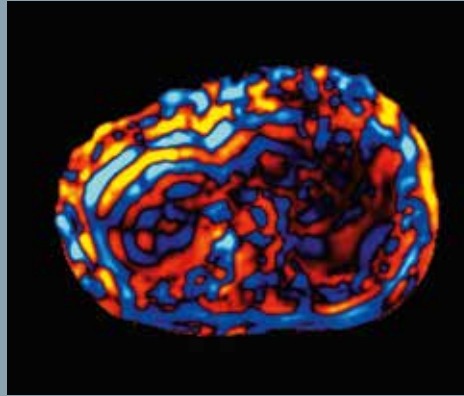

MR Elastography

MR Elastography provides the possibility to non-invasively assess variations in tissue stiffness to improve treatment decisions especially in liver fibrosis.




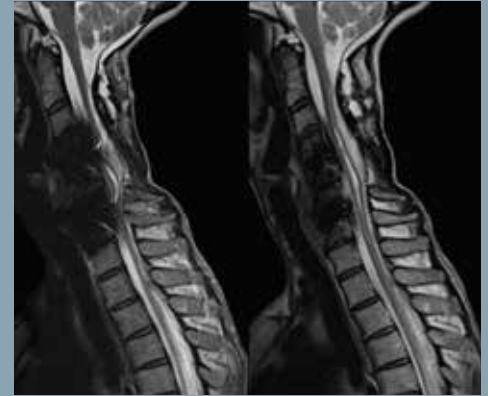


syngo WARP


syngo WARP incorporates different susceptibility artifact reduction techniques. The *syngo* MR D13 software comes with basic *syngo* WARP functionality: 2D TSE sequences are combined with high bandwidth protocols and optimized RF pulses tailored to reduce susceptibility artifacts, e.g. originating from MR conditional metal implants.

See finest anatomical detail and
visualize functional processes

Advanced clinical applications,
higher image quality

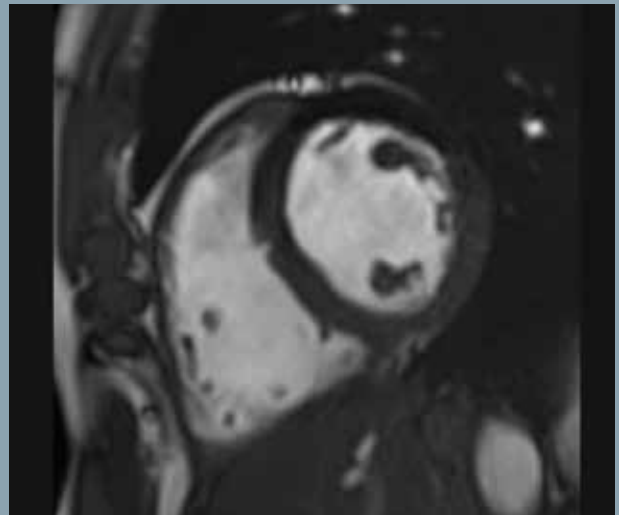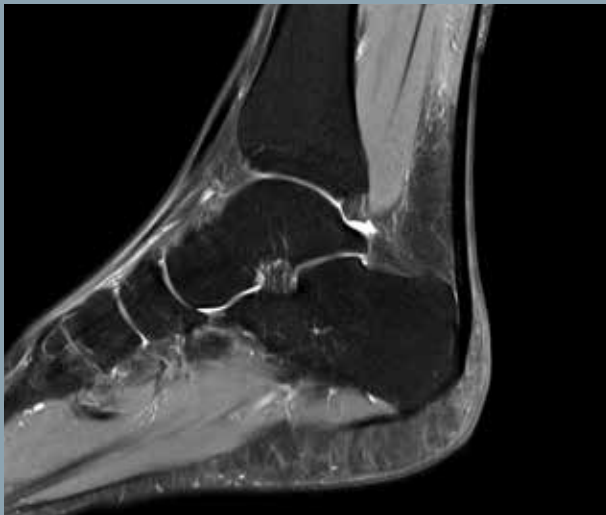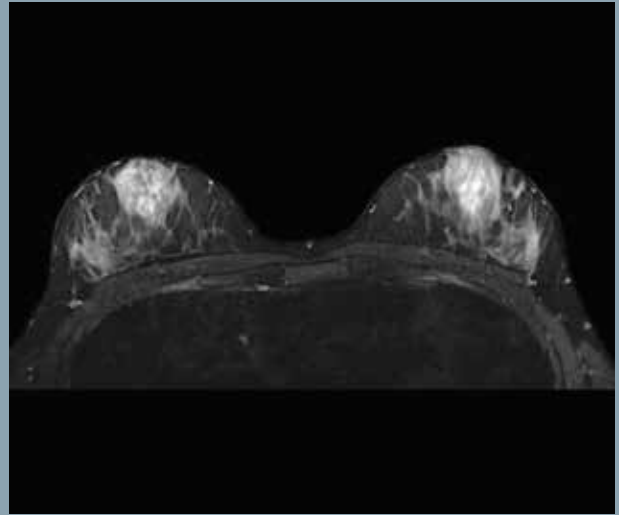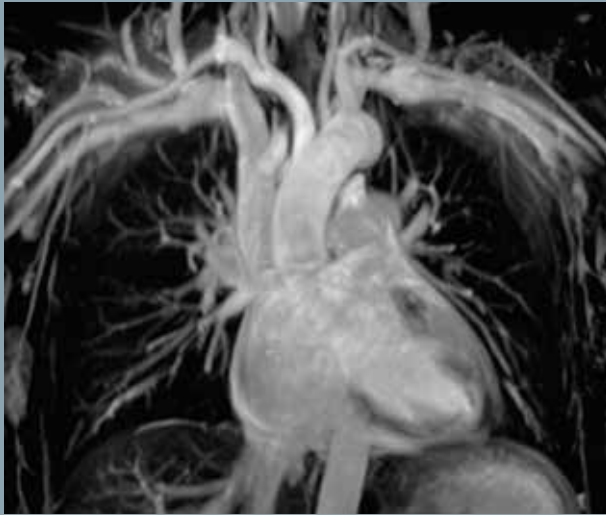


MAGNETOM Verio

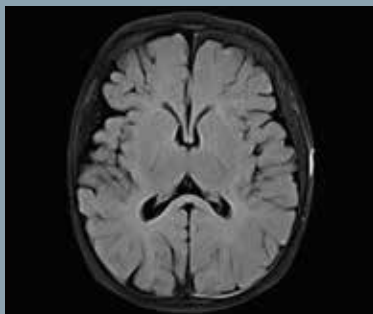

Dark Fluid, matrix 256, SL 4 mm,
TA 2:28 min

MAGNETOM Skyra^{fit}

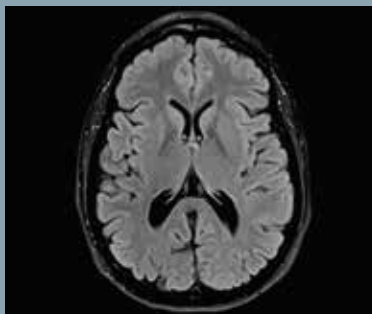

Dark Fluid, matrix 320, SL 3 mm,
TA 3:54 min

Excellent contrast

Thanks to Tim 4G, the contrast of the white and the grey matter can be improved. Consequently, more accurate diagnoses can be achieved.

*Left image:
Children's Mercy Hospital,
Kansas City, USA*

MAGNETOM Verio



PD TSE, matrix 320, TA 4:27 min

MAGNETOM Skyra^{fit}




PD TSE, matrix 384, TA 3:21 min

Shorter scan time


Due to the higher number of RF channels and the high-density coils, scan time can be reduced by up to 25% in an examination of the elbow.¹ Even if the matrix parameter is larger.

MAGNETOM Verio

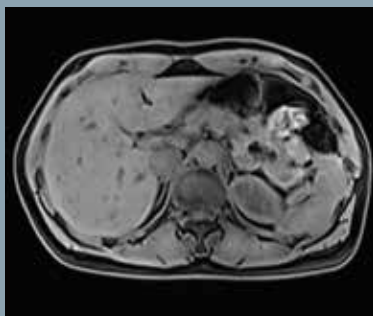

T1 3D VIBE FatSat, GRAPPA 2,
matrix 256, TA 19.12 s

MAGNETOM Skyra^{fit}

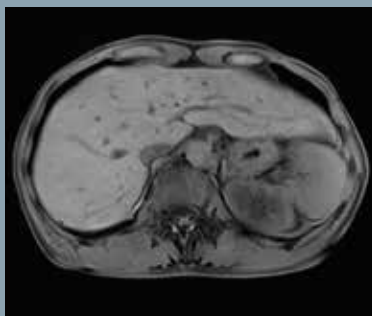

T1 3D VIBE FatSat, GRAPPA 2,
matrix 320, TA 19.38 s

Higher resolution

High-resolution imaging of the abdomen with better SNR and excellent fat saturation. In nearly the same scan time.

¹ Data on file. Results may vary.

Upgrade now

and bring your MAGNETOM Verio to the next level

How to make sure that your MAGNETOM® Verio stays with the new generation of scanners?

Upgrade it now to MAGNETOM Skyra^{fit} with Tim® 4G, the 4th generation of proven Tim (Total imaged matrix) integrated coil technology and enjoy its exceptional flexibility, accuracy, and speed.

In addition, Dot® (Day optimizing throughput), the next movement in MRI, will bring you reproducibility and consistency with a greater ease-of-use and higher productivity.

Upgrade your MAGNETOM Verio to MAGNETOM Skyra^{fit} and you'll see the extensive list of new components and applications that you will receive with the upgrade, just like having a new system.

MAGNETOM Skyra^{fit}
A Tim+Dot System

SIEMENS



New hardware:

- New RF system with 48 or 64 independent channels for faster imaging and higher signal-to-noise ratio (SNR)
- New integrated Tx/Rx Body Coil
- Fully digital with DirectRF®
- Tx/Rx real-time feedback loop for dynamic RF
- New Tim patient table for easier patient handling
- New covers

New standard coils:

The new Tim 4G coil technology with Dual-Density Signal Transfer, DirectConnect®, and SlideConnect® technology combines key imaging benefits: Excellent image quality, high patient comfort, and unmatched flexibility.

- New Head/Neck 20
- Spine 32
- Body 18
- Flex Large 4
- Flex Small 4

New optional coils available with a Tim 4G system:

- 2nd Body 18
- Peripheral Angio 36
- Foot/Ankle 16
- Hand/Wrist 16
- Tx/Rx 15-Channel Knee Coil
- Tx/Rx CP Head Coil
- 4-Channel Special-Purpose Coil

Tim 4G coils benefits:

- Designed for highest image quality combined with easy handling
- High element density of the coils that increases SNR and reduces examination times
- DirectConnect and SlideConnect technology reduce patient set-up time significantly
- Light-weight coils with an open design allow for highest patient comfort, resulting in better patient cooperation and image quality
- No coil changing with multi-exam studies saves patient set-up and table time
- All coils are time-saving “no-tune” coils

New Computer system for faster reconstruction:

- New host computer
- New Image processor

New Tim Table:

- Scan range of up to 205 cm
- Can be lowered to a minimum height of 52 cm for easier patient positioning and better accessibility for geriatric, pediatric or immobile patients
- Can be moved with two clicks into the isocenter – one click to the upmost position and one click into the isocenter

New optional Tim Dockable Table:

- Increases comfort for immobile patients, patient transport, and more
- Fast dock/undock functionality for a better patient handling
- Fits the needs for patients up to 250 kg (550 lbs)

SIEMENS



Leading.
With
MAGNETOM.

www.siemens.com/skyra-fit

Upgrade now.
MAGNETOM Skyra^{fit}

Transforming 3T productivity.

Answers for life.

Not for distribution in USA.

On account of certain regional limitations of sales rights and service availability, we cannot guarantee that all products included in this brochure are available through the Siemens sales organization worldwide. Availability and packaging may vary by country and are subject to change without prior notice. Some/All of the features and products described herein may not be available in the United States.

All devices listed herein may not be licensed according to Canadian Medical Devices Regulations. The information in this document contains general technical descriptions of specifications and options as well as standard and optional features which do not always have to be present in individual cases.

Siemens reserves the right to modify the design, packaging, specifications, and options described herein without prior notice. Please contact your local Siemens sales representative for the most current information.

Note: Any technical data contained in this document may vary within defined tolerances. Original images always lose a certain amount of detail when reproduced

Please find fitting accessories:
www.siemens.com/medical-accessories

Local Contact Information

In the USA

Siemens Medical Solutions USA, Inc.
51 Valley Stream Parkway
Malvern, PA 19355
Phone: +1 888-826-9702
Phone: +1 610-448-4500
Fax: +1 610-448-2254

In China

Siemens Medical Park, Shanghai
278, Zhouzhu Road
SIMZ, Nanhui District
Shanghai, 201318,
P.R. China
Phone: +86-21-38895000
Fax: +86-10-28895001

In Japan

Siemens Japan K.K.
Gate City Osaki West Tower
1-11-1 Osaki
Shinagawa-ku
Tokyo 141-8644
Phone: +81 3 5423 8411

In Asia

Siemens Pte Ltd
Healthcare Sector
Regional Headquarters
The Siemens Center
60 MacPherson Road,
Singapore 348615
Phone: +65 6490-6000
Fax: +65 6490-6001

Global Business Unit

Siemens AG
Healthcare Sector
Magnetic Resonance
Henkestr. 127
91052 Erlangen
Germany
Phone: +49 9131 84-0

Global Siemens Headquarters

Siemens AG
Wittelsbacherplatz 2
80333 Muenchen
Germany

Global Siemens Healthcare Headquarters

Siemens AG
Healthcare Sector
Henkestrasse 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
www.siemens.com/healthcare

Legal Manufacturer

Siemens AG
Wittelsbacherplatz 2
DE-80333 Muenchen
Germany